

# MASSAGE GUN WITH BELT



## PRODUCT INFORMATION

Gear position: 4

Interface type: DC

Force: 60-80LB

Voltage: 7.4V

Motor revolutions: 1800-3000 rpm

Battery capacity: 1500mah

Net weight of product: 0.5kg

switch: key type

shell material: ABS+painting process

Power: brushless motor 30W


Motor diameter: 35mm

Charging time: 3-4 hours

endurance: 4-6 h

gross weight of product: 0.7kg

## PRODUCT DIAGRAM

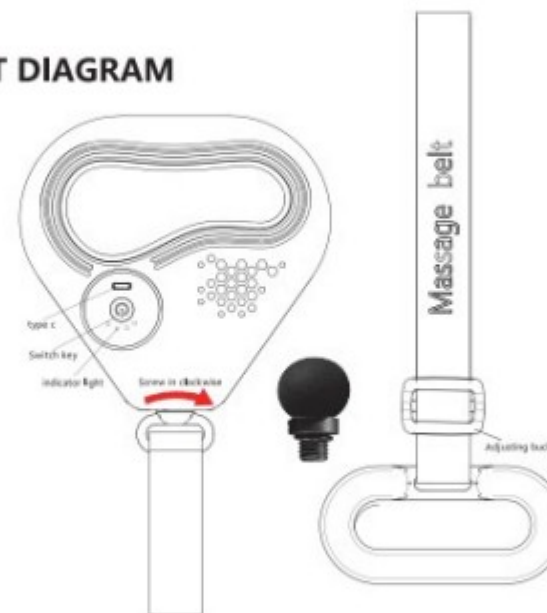

## CHARGING INSTRUCTIONS

1. Power on/off button (the indicator light is off, long press to turn on/long press to turn off)
2. Press the key gear to replace.
3. Please select 5V charging device for charging
4. For the first time, please charge for 3 hours. Please connect the data cable to the 5V adapter for charging.
5. During charging, the four lights will flicker, and the four lights will always be on when fully charged. They are red, yellow and green in turn.

## SAFETY INSTRUCTIONS:

In order to reduce the risk of electric shock, fire and personal injury, this product must be used according to the following instructions.

## FOR ADULTS ONLY

- This product can only be used on dry and clean skin surface, not through clothes. Press and gently move on the skin.
- Use the device only on soft tissues of the body without causing pain or discomfort. Do not use it on the head, any hard or bony part of the body.
- The length of massage belt is adjustable. This product does not need to be too tight when replacing the massage head.
- Bruising may occur under moderate or normal pressure. Always check the area receiving the massage.
- In case of pain or discomfort, please stop using it immediately.
- Fingers, hair or other parts of the body should be kept away from the shaft and back of the vibrating head, as they may be squeezed.
- Do not place anything in the ventilation hole of the equipment.
- Do not immerse it in water or let water enter the ventilation hole of the equipment. Do not drop the equipment or use it improperly. Only use 5 VDC to charge the equipment charger. Please carefully check the equipment and battery before each use. Do not change or modify the device in any way.
- Do not operate or charge unattended.

## MATTERS NEEDING ATTENTION

1. The interface position of this product is rubbed due to expansion and contraction. In order to prevent wear, lubricating oil should be applied. Pay attention not to get stained with clothes. The contact skin can be wiped or cleaned, but it will not be used for several times
2. This product cannot be used while charging, and you need to press and hold the ON key to use or close the key.
3. The product automatically shuts down every 10 minutes.

## The warranty permissions

The warranty period of this product is 12 months, the authority is not transferable, from the date of purchase of the user in the warranty period with valid documents can enjoy the warranty service (appearance and wear caused by mechanical operation are not within the warranty scope).

## Disclaimer :(the following is not covered by warranty)

1. Accidental damage, improper use or transportation;
2. Improper maintenance of the product;
3. Do not operate according to product instructions;
4. Unable to provide formal purchase vouchers;

## Product warranty card

Product model		Factory number		Manufacture date				
Customer data	Client name;		Zip code;					
	Detailed address;							
	Contact Number;							
Name of commodity				The cost price				
Maintenance record								

The contents of this manual and the specifications of this product are subject to change without prior notice.  
The Company reserves the right to make changes to the specifications and materials contained therein without notice, and will not be liable for damages (including consequential damages), including but not limited to typographical errors and other errors related to this publication, caused by reliance on the materials cited.  
This manual is for reference only, does not constitute any form of commitment, product configuration information, specifications and parameters are subject to the actual purchase.

## Certificate of approval

Inspection Conclusion:   
Inspector No. : \_\_\_\_\_  
Inspection Date: \_\_\_\_\_